

HOW TO USE THE G-RED (624 NM) FUNCTION

- ✓ **Sonic massage, nutrition, pore contraction, improvement of elasticity**

► STEP 1

While in standby mode, press the power button 3 times to switch to G-Red mode. The LED light will turn on and operation will start.


► STEP 2

Apply cosmetics (nourishing cream or wrinkle care products) to the skin and use according to the directions shown in the illustration.



► STEP 3

Lightly place the blade (flat side) on your skin and move it slowly in the directions shown in the illustration. Once the hand makes contact with the hand electrode, a functional LED light will turn on. The ion supply and light therapy will start.

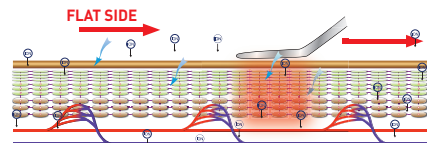

* While operating, you may feel a slight tingling due to circulating current. However, it is harmless to the human body.

* **USAGE TIME: 5 Minutes**

EFFICACY OF REJUVENA™ PREMIUM

1. WATER PEELING


When used with water, the sonic waves from Rejuvena™ Premium deeply penetrate the skin, removing impurities from the pores and surface of the skin.

2. ION

Negative and positive ions, generated from the blade, allow the penetration of cosmetic nutrients into the skin. The firmer feeling and more lifted look will make you feel remarkably younger.

3. LIGHT THERAPY

Blue light helps to improve troubled skin. Red light reduces free radicals and soothes skin damage caused by UV rays.

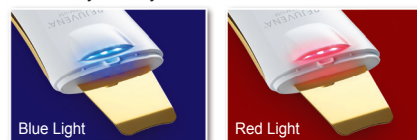

4. GOLD THERAPY

To the ancient Chinese, gold was the secret to youthful skin. Today, gold continues to be used for luxurious skin treatments.

* **Product efficacy varies from individuals.**

Model : REJUVENA™ PREMIUM
Weight : 90g
Rated voltage : DC 5V
battery : 3.7V / 450mAh (Li-Polymer)
Charging Time : 60 minutes
Size : 170 X 52 X 15.5 (mm)
Material : ABS
Country : Made in Korea
Manufacturer : RNU CORPORATION
Address : 58, Maeyeo-ro, Dong-gu, Daegu, Korea 41065
TEL : +1-855-338-6410
FAX : +1-808-597-1651

PURE C SERUM™

Our advanced Pure C Serum™ is a unique, oil-free patented Vitamin C formulation that contains powerful antioxidants designed to help reduce visible signs of aging and sagging skin.

► **Available for separate purchase**

PRECAUTIONS

- 1 Keep the Rejuvena™ Premium away from infants and young children.
- 2 Do not disassemble or modify the product.
- 3 Charge the Rejuvena™ Premium fully before the first use.
- 4 When using the ion function, never let your eyes get exposed to the functional LED for long periods of time.
- 5 Caution! If you immerse the device in water or a large amount of water gets into the device, it may cause a malfunction.

WARRANTY REGULATIONS

- 1 The Rejuvena™ Premium is subject to a 3-year limited warranty from the date of purchase.
- 2 The warranty period starts from the date of purchase, so please keep your receipt. If the date of purchase cannot be confirmed, the warranty period will be based on the date of production.
- 3 A fee will be charged if you need a repair for the following reasons. These cases are not covered by warranty.
 - A. In case of breakage, caused by violations of included user's manual instructions.
 - B. In the case of failure due to external factors such as shock and dismantling of the equipment.
 - C. In the case of failure due to natural disasters.
 - D. In case of running out of consumable parts.

RICHWAY & FUJI BIO INC.

1314 SOUTH KING STREET SUITE 520, HONOLULU, HI 96814 USA
TEL: (808) 589-2800 TOLL-FREE: (855) 338-6410
FAX: (808) 597-1651
WWW.RICHWAYANDFUJIBIO.COM

SHAOXING RICHWAY TRADING CO., LTD.

RICHWAY EUROPE B.V.
RICHWAY & LIFE JAPAN CORPORATION
R&L CORPORATION

FDA DISCLAIMER

These statements have not been evaluated by the Food and Drug Administration. This product is not intended to diagnose, treat, cure, or prevent any disease. Specific medical advice should be obtained from a licensed health care practitioner.

PUBLISHED 09092016

PEELING · ION TREATMENT · LIGHT THERAPY

REJUVENA™ PREMIUM

[User's Guide]




THE
RICHWAY
GROUP

RICHWAY & FUJI BIO INC.


1314 South King Street Suite 520, Honolulu, HI 96814 USA www.richwayandfujibio.com
Tel: 808-589-2800 Toll-Free: 855-338-6410
Fax: 808-597-1651

COMPONENTS AND PARTS

► COMPONENTS

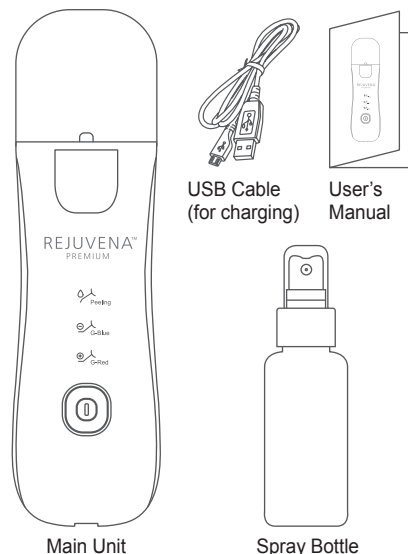


► NAMES OF PARTS AND DISPLAY SCREEN

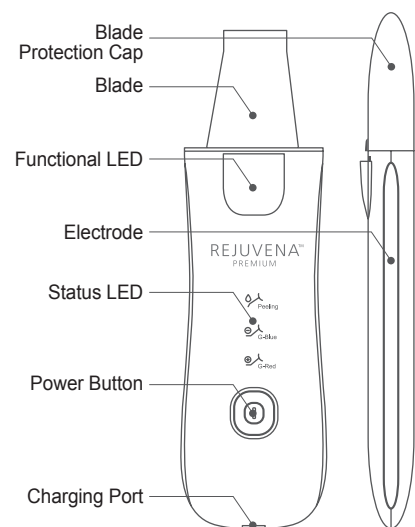

HOW TO USE THE REJUVENA™ PREMIUM

1. OPEN PROTECTIVE CAP

As shown in the picture, lightly lift the tip of the blade protection cap.

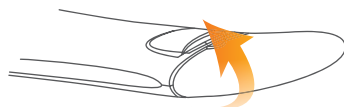

2. TURN ON THE POWER

Press and hold the power button (approximately 2 seconds). A sound will trigger and the power will turn on to standby mode. The LED lights will consecutively flash.



3. POWER OFF

Press and hold the power button (approximately 2 seconds) while operating or in standby mode. A sound will trigger and the power will turn off.

4. AUTOMATIC TURN OFF


After 5 minutes during peeling, G-Blue, or G-Red settings, the device will automatically change to standby mode. After 5 minutes in standby mode, without movement, the device will automatically turn off.

5. LOW BATTERY

If the product runs out of battery power while in use, 3 LED lights will blink 5 times and then the product will shut off automatically.

6. CHARGING

Connect the Micro USB Charger into the charging socket located on the bottom of the device. Once the device's battery begins to charge, the red LED light will turn on. The blue LED light will turn on once the product has finished charging. After the product has been fully charged, please remove the charging cable before continuing to use the product.

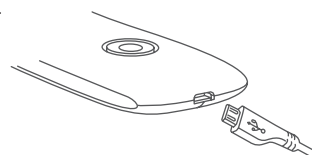

* This product cannot be used while charging.

HOW TO USE THE SOFT PEELING FUNCTION

✓ Cleaning, removal of old and dead skin cells


► STEP 1

While in standby mode, press the power button once to switch to Peeling mode. The LED light will turn on and operation will start.

► STEP 2


Cleanse the face and lightly pat with a towel. Some water should remain on the skin before use. If there isn't sufficient moisture on your skin, use after applying a light mist of water on your face. For best results, use water with a pH of 6.5.

* If there is too much moisture on the blade, it may not work effectively. Use after removing excess moisture.

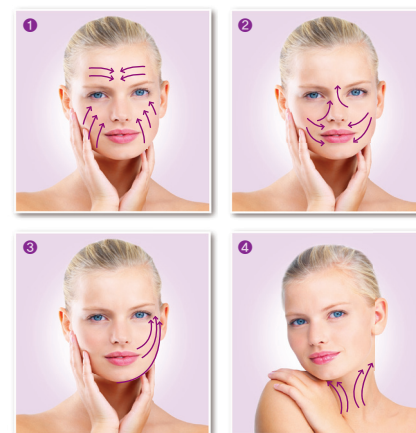

► STEP 3

Lightly place the blade (blade side) on your skin and move it slowly in the directions shown in the illustration.

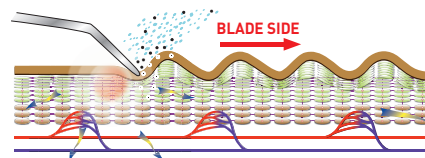

* After soft peeling, your skin may feel dry. Use a moisture essence to keep your skin hydrated.

* USAGE TIME: 5 Minutes

HOW TO USE THE G-BLUE (470 NM) FUNCTION

✓ Sonic massage, nutrition, skin trouble and deep cleaning

► STEP 1

While in standby mode, press the power button 2 times to switch to G-Blue mode. The LED light will turn on and operation will start.

► STEP 2


Apply Pure C Serum™ or your favorite cosmetic (ampoule or essence) to the skin and use according to the directions shown in the illustration.

* G-Blue mode is helpful in the absorption of vitamins.



► STEP 3

Lightly place the blade (flat side) on your skin and move it slowly in the directions shown in the illustration. Once the hand makes contact with the hand electrode, a functional LED light will turn on. The ion supply and light therapy will start.

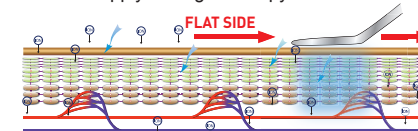

* While operating, you may feel a slight tingling due to circulating current. However, it is harmless to the human body.

* USAGE TIME: 5 Minutes